

IXplore IX85




Automated Inverted Microscope System

The IXplore™ IX85 platform delivers an unmatched level of customizability, allowing you to design or build an intelligent, high-performance imaging system that meets your specific goals. With an industry-leading field number (FN) of 26.5mm plus an array of advanced end-to-end imaging and workflow features, the IXplore™ IX85 means you can capture more than ever before while dramatically reducing your acquisition times.


Unmatched Field of View

View more of your sample in each image with a 26.5mm FN on two integrated ports. Cover larger areas in fewer images, get more datapoints from each image, minimize the need for image stitching, and speed up your workflows.

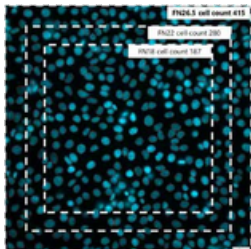


Multiple Magnifications


Our built-in magnification changer offers multiple options (1x, 1.6x, and 2x), enabling you to tailor imaging to your camera pixel size.

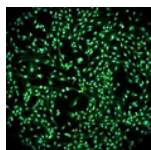

1x

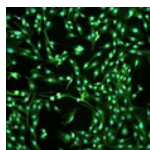


1.6x

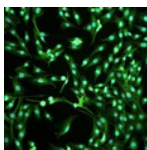

2x

New Automatic Correction Collar

Our new correction collar automatically fine-tunes your objective settings and optimizes image quality by reducing spherical aberrations caused by variations in cover glass thickness. Create sharper, more detailed images with the simple click of a button.

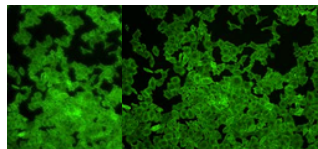

Left: without automatic correction collar adjustment.
Right: with automatic correction collar adjustment.

One Platform, Infinite Possibilities

Dual Built-In Ports

Two built-in ports (left and right) and an optional third port enable the capacity to create a multi-camera, multimodal configuration.

Open Frame Design

Add or swap modules to customize functionality. As your research needs grow you can also upgrade from one to two decks.

Built-In Threaded Holes

Threaded holes around the camera ports turn your microscope into a canvas for innovation.



IXplore IX85

Specifications

	IX85PIZF	IX85P2ZF
Microscope frame	Optical system	UIS2 optical system
	Revolving nosepiece	Motorized 6-position revolving nosepiece (DIC slider attachable) One position for Automated Correction Collar Simple water proof structure
	Focus	Stroke: 10.5 mm Minimum increment: 0.01 um Maximum nosepiece movement speed: 3mm/s
	Intermediate Magnification Changer	3 positions (Coded) 1X / 1.6X / 2X
	Light path selection	Motorized 4 positions Eyepiece 100%, left 100%, right 100%, eyepiece 50%/left 50%
	Deck insert layer	1 layer2 layers
	Maximum port field number	Left/Right side port: FN26.5 BI port: FN22 Deck right side port: FN18Left/Right side port: FN18 BI port: FN22 Deck right side port: FN18
	Dimensions (IX85 Microscope frame with illumination pillar)	540mm (D) x 321mm (W) x 667mm (H)540mm (D) x 321mm (W) x 732mm (H)
Focus compensator	TruFocus Z drift compensator	Offset method (Focus search, one-shot focus, continuous focus), Class 1 laser product, laser wavelength: 830nm
Transmitted light illuminator		Pillar tilt mechanism (30° inclination angle, with vibration reducing mechanism), Condenser holder (with with 88 mm stroke, refocusing mechanism), Field iris diaphragm adjustable, 4 filter holders Light source: High power LED light source
Observation tube	Widefield (FN22)	· U-TB190BK Wide field tilting binocular · U-B190 Wide field binocular · U-TR30-2/U-TR30NIR Wide field trinocular
Stage	Motorized stage	· IX5-SSA: Stage stroke: X: 116mm x Y: 78mm, maximum stage movement speed: 40mm/s. Knob controller · 3rd party motorized stage
	Mechanical stage with left handle Mechanical stage with right handle	Stage stroke: X: 116mm x Y: 78mm, stage position locking function
	Right handle stage	Stage stroke: X: 50mm x Y: 50mm
	Gliding stage	Upper circular stage 360° rotatable, 20 mm (X/Y) travel
	Plain stage	232 mm (X) x 240 mm (Y) stage size, stage insert plate exchangeable (ø110 mm)
Condenser	Motorized long working distance condenser	W.D. 27 mm, NA 0.55, motorized turret with 7 position slots for optical devices (3 positions for ø30 mm and 4 positions for ø38 mm), motorized aperture and polarizer
	Long working distance universal condenser	NA 0.55, W.D. 27 mm 5 positions for optical devices (3 positions for ø30 mm and 2 position for ø38 mm)
	Ultra long working distance	NA 0.3, W.D. 73.3 mm, 4 positions for optical devices (for ø29 mm)
Fluorescence illuminator	L-shape-fluorescence illuminator	L-shaped design with exchangeable FS and AS modules, slider shutter and ND filter pocket
Fluorescence mirror turret	Motorized fluorescence mirror turret	Motorized turret with 8 positions, built-in shutter, simple waterproof structure
Fluorescence light source		· U-LGPS: LED and LDP light source, Class 1 laser product · 3rd party LED light source
Control unit (IX5-MCZ)		Nosepiece position, light path selection, filter turret position, FL shutter ON/OFF, DIA LED power, DIA LED ON/OFF, 4 customizable buttons
Control box	PC interface	USB (Type-C), RS-232C
Operating environment	- Indoor use - Ambient temperature: 5o to 40o C (41o to 104o F) - Maximum relative humidity: 80% for temperatures up to 31 oC (88 oF), decreasing linearly through 70% at 34 oC (93 oF), 60% at 37o C (99o F), to 50% relative humidity at 40o C (104o F) - Supply voltage fluctuations: Not to exceed ±10% of the normal voltage¹⁾	



EVIDENT
OLYMPUS

LABEX
FLUOCCA
IMARIS

GREINER
BIOSCIENCE

APPLIED
SPECTRAL
IMAGING

SAKURA
FINETEK

ZTEC
CRAFTK

TEKSQRAY
FEATHER

DYNAMICA
FROILABO

FASTER
ARCTIKO
KUHNERT
HETTICH

HBK
CAMECA
IX CAMERAS

THERMO EM
CRYSTAL

IX Applications for ICSI/IMSI

Optimized for a Smoother ICSI Workflow

IX73 Microscope

Motorized Observation Method Switching

- Quickly change observation methods by pressing a button on the hand switch
- Spindle observation is standard with a semi-motorized system; IMSI units are available




IX83 Microscope

Fully Motorized and Highly Advanced

- Standard configuration supports IMSI and spindle observation
- Quickly change observation methods and magnification using the controller
- Focus position and brightness automatically adjust when you change magnifications
- Easily operate the focus and switch between camera and visual observation using the remote control, reducing sample vibration



Applied Technologies

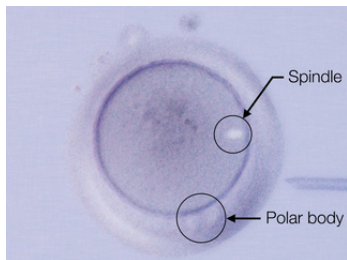


Easy Oocyte Condition Checks through Spindle Visualization

Olympus' observation method enables you to easily view the metaphase II (MII) spindle through the microscope oculars



Relief Contrast

Olympus' relief contrast enables 3D imaging of oocytes so you can check the condition of the zona pellucida

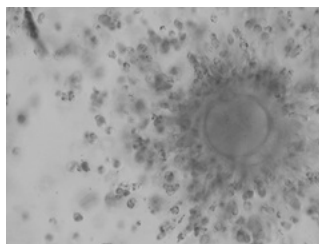

Polarization

Observe the birefringence of spindles using polarized light