


CKX53

Cell Culture Microscope

The CKX53 inverted microscope delivers stable performance and a comfortable workflow for various cell culture needs, including live cell observation, cell sampling and handling, image capture, and fluorescence observation. The system provides clear, reproducible, and high contrast images with a wide visual field.

Specifications

Illumination		Brightfield	Phase Contrast Entry	Phase Contrast Standard	Fluorescence
Focus		Revolving nosepiece vertical movement system using the coarse and fine focusing knobs. Stroke: 20 mm (Focal point: up to 18.5 mm from the plain stage top surface) Stroke per rotation : 36.8 mm (Coarse), 0.3 mm (Fine)			
Stage	Plain stage	252 mm (D) × 200 mm (W) Exchangeable transparent insert plate is incorporated			
	Mechanical stage	Options		XY coaxial knob place on right stageside of the plain stage Microplate holder equipped with the escape function stage stroke: X = 110 mm, Y = 74 mm	
	180 mm (D) × 70 mm (W)				
	Substage				
Condenser		Maximum numerical aperture: 0.3 Working distance: 72 mm Applicable objective magnification 2X, 4X, 10X, 20X and 40X up to 190 mm height tissue flask can be loaded on the stage without detachable condenser			
Observation tube		Fixed Trinocular tube, inclined 45 degrees Interpupillary distance 48 – 75 mm Light path: eyepiece/camera port = 100/0 ⇄ 0/100			
Dimensions (W × D × H)		200 mm × 498 mm × 454 mm (Phase contrast entry configuration)			
Weight		6.9 kg	7.1 kg	8.4 kg	12.5 kg

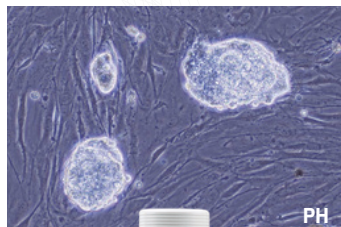

CM30

Incubation Monitoring System

The CM30 incubation monitoring system provides quantitative data remotely— the system will periodically scan it, count the number of cells, and determine confluency. So you can monitor your cultures' progress without entering the clean room.

Specifications

Installation environment (inside the incubator)	Temperature: 37 °C (98.6 °F) ± 0.3 °C (0.5 °F), humidity: 0—99%
Applicable vessels*	Petri dish (90 mm (3.54 in.), 100 mm (3.94 in.)) Microplate (6 well, 12 well, 24 well, 48 well, 96 well) Flask (T25, T75, T80, T150, T175, T225) Multi-layer flask
Optical performance	Field of view (H × V): 2.84 mm × 2.13 mm (0.11 in. × 0.08 in.); (image size per one shooting) Image size: 1280 × 960 pixels Illumination wavelength: λ = 630 nm (LED) Illumination method: epi-oblique illumination
Cable length	Approx. 4.5 m (14.8 ft)
Sterilization resistance	Autoclave sterilization (for vessel holder and sponge rubber only) UV ray sterilization
Disinfection	Hydrogen peroxide (H ₂ O ₂) gas sterilization (CM30H only) Peracetic acid disinfection (cold sterilant)
Weight	Approx. 3.1 kg (6.8 lb)

IPSC (Induced pluripotent stem cells)



See more in our Application Gallery

EVIDENT
OLYMPUS

LABEX
FLUOCCA
IMARIS

GREINER
BIOSCIENCE

APPLIED
SPECTRAL
IMAGING

SAKURA
FINETEK

ZTEC
CRAFTK

TEKSQRAY
FEATHER

DYNAMICA
FROILABO

FASTER
ARCTIKO
KUHNERT
HEITICH

HBK
CAMERA
IX CAMERAS

THERMO EM
CRYSTAL