




EVIDENT
OLYMPUSLABEX
FLUOCCA
IMARISGREINER
BIO SCIENCEAPPLIED
SPECTRAL
IMAGINGSAKURA
FINETEKZTEC
CRAFTTEKTEKSQRAY
FEATHERDYNAMICA
FROLLABOFASTER
ARCHIKO
KUNHIERHBK
CAMERA
IX CAMERASTHERMO EM
CRYSTAL

SZX7

Stereo Microscope System


The SZX7 stereo microscope's Galilean optical system offers excellent image quality, especially when using a digital microscope camera. The microscope is designed for life science imaging applications with high color fidelity optics, a 7:1 zoom ratio, and a universal LED stand.

Observation Images

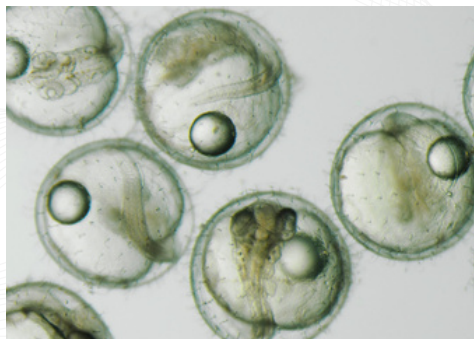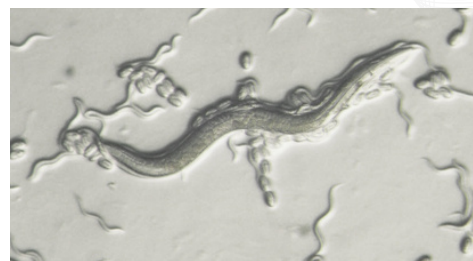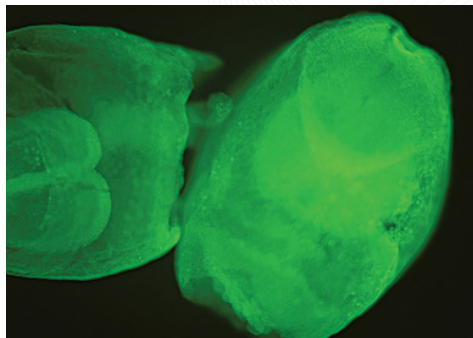

Specifications

Optical system	Galilean type optical system
Zoom microscope body	Zoom range 0.8X – 5.6X (zoom ratio 1:7), Lead-free
Observation tube	1. SZX-BI45: 45° binocular, 45° inclination 2. SZX2-TR30: 30° trinocular, 30° inclination 3. SZX2-TR30PT: 30° trinocular, 30° inclination 4. SZX2-LTTR: Ergonomic long tilting trinocular, 5 – 45° variable inclination All observation tubes: Lead-free Interpupillary distance adjustable range: 50 – 76 mm
Extendable eyepoint adjuster	SZX2-EEPA: Height adjustment range: 30 – 150 mm (with a scale attached)
Objectives	DFPL0.5X-4, DFPL0.75X-4, DFPLAPOIX-4, SZX-ACHIX, DFPLAPOI.25X, SZX-ACHI.25X-2, DFPL1.5X-4, DFPL2X-4
Eyepiece	WHSZ series
Dimensions (W × D × H)	194 mm × 253 mm × 403 mm (Standard Set Configuration)

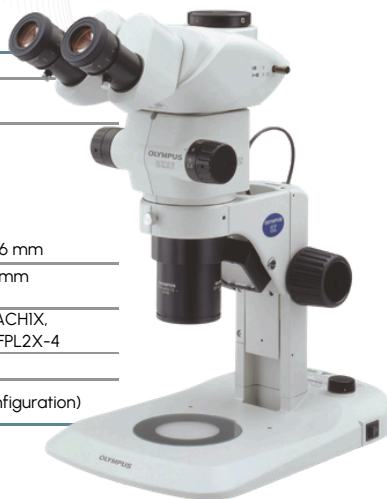

SZX/SZ Accessories



SZ2-ILST

LED illuminator stand

This LED stand features a thin design to keep sample positions low and optimize operability. Simultaneous transmitted and reflected light are available on this stand. LED light offers a long lifespan and constant color temperature at any intensity.



SZ2-ILA

Transmitted illumination attachment

Used with the SZ2-ST stand, this cost-effective illumination stand offers bright, uniform illumination from low to high magnifications. A tiltable mirror provides direct and oblique illumination for low-contrast specimen. The available 6 W and 100 W lamphouses deliver the necessary power for various illumination needs.



SZ2-CLGDI

Dual interlock light guide

Standard oblique semi-rigid fiber optic light guide. The light source position on the rear side of the stand saves desk space.



SZ2-CLGSF

Flexible light guide

A single fiber optic guide is fixed at the back of the objective so as not to disturb microscope operation.



SZX2-ILLTQ/SZX2-ILLTS

LED Transmitted Light Illumination bases

SZX2-ILLTQ with quad-position turret enables you to choose cartridges and to switch between brightfield (standard/high/low), oblique (standard/high/low), darkfield, polarized illumination, and the shutter with a simple turn. A one-position LED illumination base is also an option (SZX2-ILLTS)



SZX-SDO2

Side-by-Side Discussion Tube

Ample distance (650 mm) between the main and secondary observer for easy imaging without disturbing the microscope operation. The color of the built-in pointer can be selected to contrast the specimen.

EVIDENT
OLYMPUS

LABEX
FLUOCCA
IMARIS

GREINER
BIO SCIENCE

APPLIED
SPECTRAL
IMAGING

SAKURA
FINETEK

ZTEC
CRAFTTEK

TEKSQRAY
FEATHER

DYNAMICA
FROLABO

FASTER
ARCTIKO
KÜHNER
HETTICH

HBK
CAMERA
IX CAMERAS

THERMO
EM
CRYSTAL